


Installation instructions

To run, Holter LX software must be installed on your hard disk. No other Windows programs should be running when you install the software.

Note: If you already have Holter LX software installed and want to update to the newest version, 5.2, you can start at Step 4. You must install the new Java version, 1.4.2, for Holter LX version 5.2 to run.

To install the Holter LX software:


1. Insert the NorthEast Holter LX Software CD into your computer's CD drive. A window opens on your monitor acknowledging the CD and displaying icons indicating what is present on the CD. They include one for Adobe Acrobat Reader (its label starts "Adb"), one for Java 2 Runtime Environment (its label starts "j2re"), one for the Rainbow key that is specific to your system (its label starts "rainbow"), and one for NorthEast Holter LX (its label starts "nem"). In addition, there are document files containing a copy of this and other NorthEast product operator's manuals, and a "readme" text file.

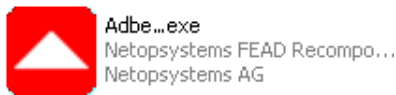

Acrobat Reader icon

(We recommend that you print and read the readme file as soon as installation is complete.)

2. Double click on the Adobe Acrobat Reader icon and follow the installation instructions as they appear on the screen, pressing the Enter key whenever prompted to go to the Next window. When completed, the CD window re-appears.

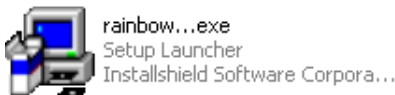


Rainbow icon

Enter key whenever prompted to go to the Next window. When completed, the CD window re-appears. Plug the software key back into its port.

Note: If you are asked to reboot after Rainbow installation, please reboot and then continue with step 4.

4. Double click on the Java icon and follow the instal-

lation instructions as they appear on the screen, pressing the Enter key whenever prompted to go to the Next window. When completed, the CD window re-appears.

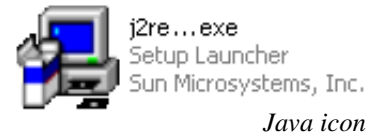

Note: If you are asked to reboot after Java installation, please reboot and continue with step 5.

5. Double click on the NorthEast icon and follow the installations instructions as they appear on the screen. Each time you are prompted to go to the Next instruction window, press the Enter key to do so. When you are prompted to select Finish, press the Enter key to do so. After a short delay, the Setup window opens.



NorthEast Monitoring icon

6. Type appropriate information in each of the fields, including your name and the name of your facility. In the five fields for the Report Heading, type the five lines you want to appear at the top of the front page of each Holter report. Leave the Number of patients set to 30. Click on the radio button next to "Digital recorder with a software key," then carefully type your registration number in the Registration field. When the information is complete, click OK to close the Setup window.

Setup window

7. Close the CD window. The Holter LX software launches and is ready to be used. To exit the LX software, select Patient > Exit or click the Close button.